

OPERATING INSTRUCTIONS

POWER

PWR BUTTON: Press the PWR button to turn the unit On. To turn Off: press the PWR button again. The control wand is programmed to remember its last function prior to being turned off. Zone 1 in SELECT mode, will be activated first if the unit was unplugged.

PWR button Timer: (15 or 30 min) When the unit is turned On the 15 minute massage Timer will activate and a green light will appear. To activate the 30 minute Timer, press and hold down the PWR button until the amber light appears.

Other ways to "turn on" PWR: Press any Zone to activate the system.

HEAT: The HEAT function may be used with massage or Independently of the massaging functions. To activate simply press the LO button for low heat or the HI button for high heat.

A separate 30 minute timer is activated when heat is turned on.

ZONES / MODES

WAVE: Press the WAVE button and the massage transducers will travel down and back up through the programmed Zones (1-3), creating a wave-like effect.

To activate: Press the WAVE button.

To program: Press any or all Zones (1-3) to activate the desired areas of massage for the wave-like action. To deactivate any zone, simply press the associated button.

PULSE: Press the PULSE button and the Zones selected will simultaneously pulse on and off.

To activate: Press the PULSE button.

To Program: Press any or all Zones(1-3) to activate the desired areas of massage for the pulse action. To deactivate any zone, simply press the associated button.

SELECT: Press the SELECT button to allow any or all zones to operate continuously.


To activate: Press the SELECT button.

To program: Press any or all Zones(1-3) to activate the desired areas of massage. To deactivate any zone, simply press the associated button.

SPEED / INTENSITY

SPEED: Controls the rate at which the massage action travels. The Speed feature operates in the WAVE or PULSE modes. Quick presses will generate a slight change. Holding down either the + (Up) or - (Down) button will achieve a greater increase or decrease.

INTENSITY: Controls the Intensity of the motors while in any function. Quick presses will generate a slight change. Holding down either the + (Up) or - (Down) button will achieve a greater increase or decrease.

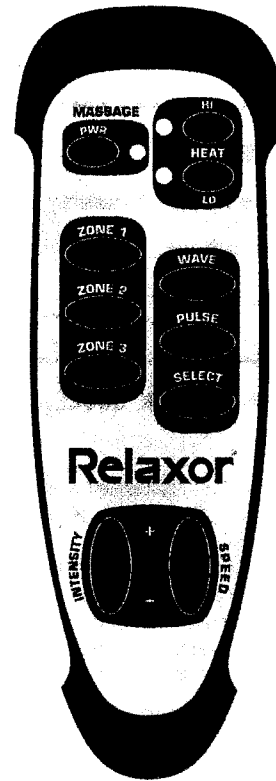

Relaxor[®] Massage System

OWNERS MANUAL and OPERATING INSTRUCTION

Dear Customer:

Welcome! You are now the proud owner of a state-of-the-art **Relaxor** massage product.

Your massage system uses a proprietary microprocessor that selects unique massage patterns and heat from a memory bank menu to provide a soothing massage to tired, knotted muscles, and relieve the tensions of stressful living with an overall dimension of quality relaxation.

This unique massage product is a scientifically designed electronic, electro-mechanical system with multiple pulsing transducers, strategically located to provide a travelling wave-like massage over your body at selectable speeds and intensities.

You now have total relaxation at your command through finger-like pressure and with soothing heat, you have the ultimate winning combination for your comfort.

We know you will enjoy your new **Relaxor** massage product.

IMPORTANT SAFETY INSTRUCTIONS

Read all instructions before using the **Relaxor** Massage System.

DANGER!

To reduce risk of electric shock.

1. Always unplug the product from the electrical outlet immediately after using and before cleaning.
2. Do not reach for a product that has fallen into water. Unplug immediately.
3. Do not use while bathing or in a shower.

WARNING!

To reduce risk of burns, fire, electric shock, or injury to persons.

OPERATE YOUR SYSTEM ACCORDING TO THIS MANUAL

USE ONLY THE POWER DEVICE SUPPLIED WITH THIS PRODUCT

1. Your **Relaxor** product should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate the **Relaxor** product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact your dealer for instructions on returning it for examination and repair.
4. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
5. Do not use outdoors.
6. Keep the cord away from heated surfaces.
7. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
8. Use heated surfaces carefully. May cause serious burns. Do not use heated surfaces over insensitive skin areas or in the presence of poor circulation.

MAINTENANCE INSTRUCTIONS:

No user maintenance of this product is required. If the System does not function, be sure...

- (a) The TRANSFORMER is plugged into an electric outlet.
- (b) The TRANSFORMER cord and COMMAND WAND cord are plugged into, and fully engaged, in their respective receptacles.

If there is malfunction of any kind, contact your dealer for instructions on returning it for examination and repair.

SAVE THESE INSTRUCTIONS!

© COPYRIGHT 1997 JB RESEARCH, INC.

PWR POWER

- ▲ Turn On: Press PWR
- ▲ Turn Off: Press PWR again

15 and 30 minute Timers:

- ▲ Green light = 15 minute Timer
- ▲ Hold down PWR: Amber light = 30 minute Timer
- ▲ To restart after Timer shuts off: Press PWR
- ▲ Unit may also be turned on by pressing any ZONE (1-3) key

ZONES

1	Upper Back
2	Lower Back
3	Hip Region

HEAT

Press LO or HI to turn On

HI ▲ Works with massage or independently of massage

LO ▲ 30 minute Timer

▲ To continue after Timer shuts off, Press LO or HI

3 MODE KEYS

WAVE

Massaging motors travel down and up in a wave-like action

- ▲ Press WAVE
- ▲ Select any or all Zones(1-3)
- ▲ Press Zone key (1-3) to cancel

PULSE

Massaging motors pulse On and Off

- ▲ Press PULSE
- ▲ Select any or all Zones(1-3)
- ▲ Press Zone key (1-3) to cancel

SELECT

Selective massage area

- ▲ Press SELECT
- ▲ Choose any or all Zones(1-3)
- ▲ Press Zone key (1-3) to cancel

INTENSITY

+ Controls INTENSITY during all functions

- Press +(Up) or -(Down) key for desired effect

SPEED

+ Controls traveling SPEED of massaging motors in WAVE and pulsing rate in PULSE

- Press +(Up) or -(Down) key for desired effect